

1,728 Skin Care SLDs, Now enjoy in Self-mode!



LED therapy by enjoying Self mode

In the skin clinic, as manual mode,
in the home, as self-mode



All-in-one Care

All-in-one Care from Facial Care to Body Care




Safe skin care for pain and side effects

Safe treatment through LLLT(Low Level Laser Therapy)



'SLD' is a differentiating LED technology

Reducing the angle of light irradiation
by up to 15° for stronger light penetration

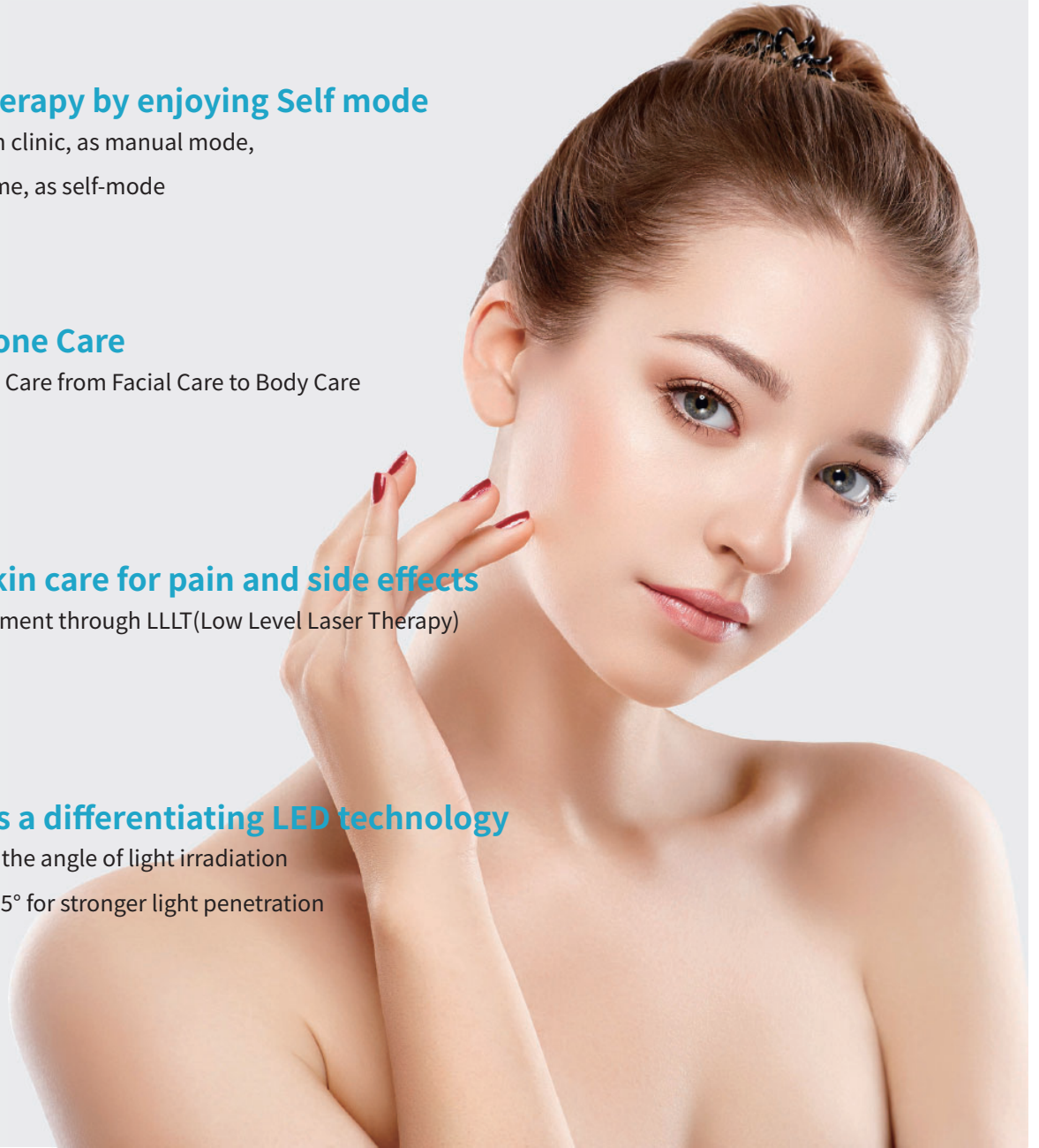

- **Sterilization** Keratinocytes, Acne
- **Circulation of Blood** Bruise, Edema, Vascular, Post-Plastic Surgery Treatment
- **Cell Reproduction** Skin Recovering, Wound Healing
- **IR Reduction of Inflammation** Anti-inflammation, Pain Relief



Acne care without pain or side effects

Blue(420nm) wavelength, which is good for sterilization can care acne effectively without treatment of any injection or medicine that can make pain or side effects.



Soothing skin after laser treatment

Red(635nm) wavelength, which is good for skin rejuvenation, helps to soothing stimulated skin and reducing flush, and IR(830nm) reduces pain and inflammation after treatment.



Anti-aging starts with LED light

Red(635nm) wavelength improves skin more tightly by making skin layer healthy, helps dull skin tone more brightly, and makes younger skin by collagen regeneration.

Before & After



Before



After



Before



After

LED Therapy with ESTHELUX

◆ Recommended distance: 10~15cm from the face and 10cm from the body.



Face(Basic)

Wavelength : Red + IR
Power : Lv.5
Time : 20 min a day



Face(Trouble Care)

Wavelength : Blue
Power : Lv.3
Time : 5 min a day



Body

Wavelength : Red + IR
Power : Lv.5
Time : 20 min a day



Scalp & Hair

Wavelength : Red + IR
Power : Lv.5
Time : 20 min a day



Face & Body

Wavelength : Yellow
Power : Lv.5
Time : 15~20 min a day

